truelife SonicBrush Junior

EN – USER MANUAL


Pack Contents

- TrueLife SonicBrush Junior
- 2x Junior heads
- Travel case white

- Charging Base
- Colored stickers
- Manual and Safety Instruction

Product Layout

1. Removable sonic toothbrush head
2. Sonic toothbrush body
3. On / Off button and changing mode
4. Mode indicators
5. Charging indicator
6. Charging base
7. Travel case

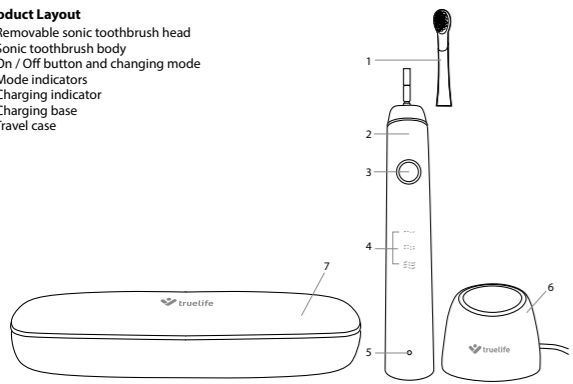Instructions for Use

1. Place the cleaning head on the brush body.
2. Moisten the head and apply toothpaste to it.
3. Place the toothbrush on the teeth and then turn it on.
4. Press the button to select the cleaning mode. When doing this, press the button for less than 5 seconds or the brush will turn off. Choose one of three modes:
 - **SLOW** – This mode is suitable if we want to teach the child how to use the sonic toothbrush correctly.
 - **MEDIUM** – This mode is for younger children who already know how to use the sonic toothbrush.
 - **FAST** – The mode with maximum power is designed for older children who already brush their teeth completely on their own.
5. Move the bristles back and forth slowly and slowly move between teeth. At the same time, gently push the brush against the tooth.
6. For maximum efficiency, the cleaning is divided into 4 sections. Each section is given 30 seconds. After a 30 second cycle, the brush pauses and alerts you to move to the next section. After 2 minutes the brush switches itself off.
7. Rinse the head and body of the brush under running water and then wipe them with a cloth or a towel.
8. Place the body back into the base or store it in a travel case.

Charging

The TrueLife SonicBrush Junior is waterproof and safe to use and recharge in the bathroom. Never, however, immerse it under water.

Plug into an electrical socket and insert the sonic toothbrush into the charging base. The charging indicator lights up. If you are charging your device for the first time, full charging will take 16 hours. Further charging of the discharged battery to fully charged will take 12 hours. In normal use, a fully-charged battery lasts 30 days. The red light on the battery indicator lights up when the battery is running low. Whilst charging, the indicator blinks red and turns green when it is fully charged.

To save energy, the charging base can be disconnected from the power supply or left on to charge the device regularly. The battery is at no risk of damage if the sonic toothbrush remains in the charging base. If the device is not used for a long time, it is normal that the indicator does not light up immediately. The indicator lights up within 30 minutes.

Maintenance

Before cleaning, make sure the appliance is not plugged in.

- Base:
Wipe the dirt from the charging stand with a soft cloth. Do not immerse the charging base in water.

- Brush:
1. Remove the cleaning head.
2. Use a soft cloth to wipe the dirt off all surfaces and return the brush back to the charging base. Do not immerse the entire body of the brush in water.
3. Rinse the cleaning head and bristles after each use.
4. Do not put the brush or its parts in the dishwasher.

- Travel case:
Wipe the surface and inside of the travel case with a damp cloth.

Recommendation: For hygiene reasons, we recommend changing the brush head at least every 3–4 months.

Specifications

Material	ABET ABS/PP/DuPont
Voltage input	DCSV 1A
Battery capacity	800 mAh
Battery life	28–30 days (used twice a day for 2 minutes)
No. of movements	56.000/min
Dimensions	26x26x231 mm
Weight	175 g

Manufacturer's Responsibility

TrueLife products are subject to a warranty for manufacturing defects. It is the responsibility of elem6 s.r.o. to repair or replace the defective part or product for a similar part or product at its sole discretion in the event of such defects. Conversion, alteration or other unauthorized tampering of the product exceeds the scope of this warranty. This is even the case if a fault, negligence, or other error in the product results in its loss or damage. The warranty does not apply to damage resulting from use, accident or common wear. Elem6 s.r.o. assumes no responsibility for any accidents, injuries, death, loss or other requirements related to the product or resulting from the use of this product. Under no circumstances is elem6 s.r.o. liable for incidental or consequential damages in connection with or arising out of the use of this product or any of its components. Repayment or replacement of spare parts and/or products may be subject to charges for transportation, handling, replacement and/or replenishment.

Print errors and information changes reserved.

SAFETY WHEN WORKING WITH THE DEVICE:

1. Carefully read the safety information leaflet.
2. For the correct handling of the product, read the current version of the manual, which is available for download at www.truelife.eu.
3. Print errors and information changes reserved.

Manufacturer:
elem6 s.r.o., Braškovská 15, 16100 Prague 6,
Czech Republic

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