 **wellneo**[®]

SONIC TOOTHBRUSH

HEALTHY TEETH AND GUMS

High frequency sonic vibrations

MODEL NUMBER: TB8061


User Manual

EN.....	3
AL	5
BG.....	7
BiH.....	11
CZ.....	14
EE	16
HR.....	18
HU	22
KS.....	25
KZ.....	27
LT	30


LV.....	32
MD.....	35
ME.....	37
MK	29
PL	42
RO.....	45
RS.....	47
RU	49
SI	52
SK	55
UA	57

(EN) PICTURES / (AL) FOTOGRAFIN / (BG) ИЛЮСТРАЦИИ / (BiH) SLIKE / (CZ) OBRÁZKY / (EE) JOONISED / (HR) SLIKE / (HU) KÉPEK / (KS) FOTOGRAFIN / (KZ) СУРЕТТЕР / (LT) PAVEIKSLIAUKAI / (LV) ATTĒLI / (MD) IMAGINI / (ME) SLIKE / (MK) СЛИКИ / (PL) RYSUNKI / (RO) IMAGINI / (RS) SLIKE / (RU) RYSUNKI / (SI) SLIKE / (SK) OBRÁZKY / (UA) МАЛЮНКИ


PICTURE 1

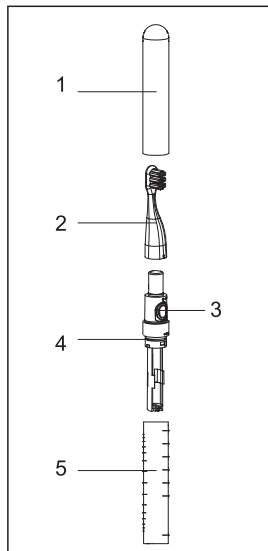PICTURE 2

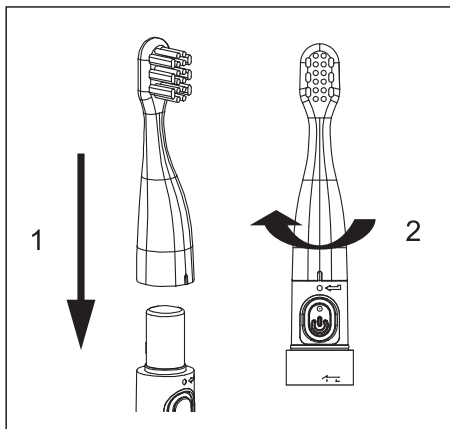PICTURE 3

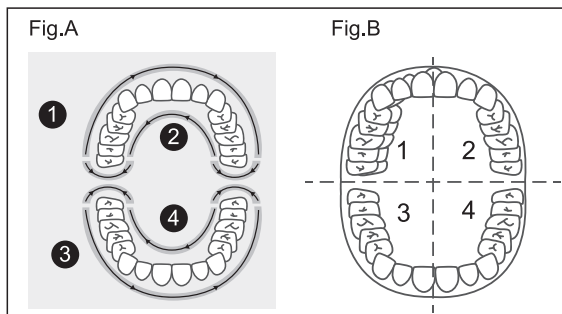Read this manual thoroughly before using and save it for future reference.

WARNINGS

1. This device is not a toy and should not be used by children or stored accessible to them.
2. This device is intended for household use only.
3. Only the batteries with a specification of 1.5V AAA can be used.
4. The instructions marked on the battery must be followed.
5. When replacing the battery, turn off the appliance first, then wipe out the water on the body. The battery must be installed correctly (polarity) and firmly.
6. Exhausted batteries are to be removed from the appliance and safely disposed of.
7. If the device is to be stored unused for a longer period (i.e. two weeks), the battery should be removed.
8. Do not use the battery which is out of date.
9. If the device is not working as it should, has received a sharp blow, has been dropped, damaged, left outdoors, or dropped into water, do not use it.
10. Do not attempt any modification or repair by yourself and ensure that any repair is conducted only by appropriately qualified technician.
11. In case of malfunction and/or damage do not attempt any repair by yourself and ensure that the product is properly disposed in accordance with local laws.
12. Do not immerse the toothbrush into the high-temperature water.
13. This device can be used by children aged from 8 years and above as well as by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the device in a safe way and understand the hazards involved. Children should not play with the device. Cleaning and user maintenance should not be done by children unless they are older than 8 and supervised.
14. Do not press the tooth or gum by the toothbrush with excessive force and avoid brushing the same position for long time to avoid any injury.
15. Please consult the doctor whether you can use the toothbrush.
16. Oral bleeding may occur when using the toothbrush for the first time or improperly. If it lasts for more than a week, please consult the doctor for suggestion.

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