

BRAUN

No touch + forehead


Thermometer



NTF 3000

GB	English 2
CZ	Česky 11
HR	Hrvatski20
HU	Magyar29
PL	Polski38
RO	Română47

RU	Русский56
SA	اللغة العربية 73
SI	Slovenščina 74
SK	Slovenčina 83
TR	Türkçe 92

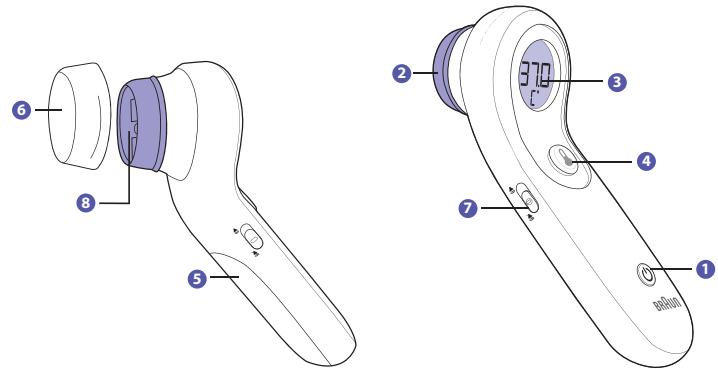

This product is manufactured by Kaz Europe Sàrl under a license to the 'Braun' trademark. 'Braun' is a registered trademark of Braun GmbH, Kronberg, Germany.

 Kaz Europe Sàrl
Place Chauderon 18
CH-1003 Lausanne
Switzerland



Printed in China

© 2015, Kaz Europe Sàrl

www.hot-europe.com

Patents: www.kaz.com/patents/braun

Made in China

NTF3000EE
P/N: 311MNTFE191R0
REV. 20NOV15

Thank you for purchasing the Braun No touch + forehead thermometer (NTF 3000). This thermometer is a high-quality product incorporating the latest technology and tested in accordance with international standards. With its unique technology, the product can provide a stable, accurate reading with each measurement. The thermometer performs a self-test every time it is switched on to always guarantee the accuracy of measurements. The No touch + forehead thermometer is intended for the intermittent monitoring of human body temperature in the home. It is intended for use on people of all ages except pre-term babies or very small (small-for-gestational age) babies.

Please read these instructions carefully before using this product and keep the instructions and the thermometer in a safe place.



WARNINGS AND PRECAUTIONS

Keep out of reach of children under 12 years.

Never use the thermometer for purposes other than those it has been intended for. Please follow the general safety precautions when using on children.

Never immerse the thermometer into water or other liquids (not waterproof). For cleaning and disinfecting please follow the instructions in the "Care and cleaning" section.

Do not store this thermometer in temperature extremes below -25 °C or over 60 °C (below -13 °F or over 140 °F) or in excessive humidity (above 95% non-condensing relative humidity).

If thermometer is stored in a location that is cooler or warmer than where it is being used, let it sit in the patient's room for 10 minutes before taking the measurement.

Do not use the thermometer if there are signs of damage on the scanner or on the thermometer itself. If damaged, do not attempt to repair the product.

Never insert a sharp object into the scanner area or any other open surface on the thermometer.

This thermometer consists of high-quality precision parts. Do not drop the instrument. Protect it from severe impact and shock. Do not twist the instrument or the measuring sensor.

This thermometer is intended for household use only.

Use of this thermometer is not intended as a substitute for consultation with your physician.

Temperature elevation may signal a serious illness, especially in neonates and infants, or in adults who are old, frail, or have a weakened immune system. Please seek professional advice immediately when there is a temperature elevation and if you are taking temperature on:

- Neonates and infants under 3 months (Consult your physician immediately if the temperature exceeds 37.4 °C [99.4 °F])
- Patients over 60 years of age (Fever may be blunted or absent in older patients)
- Patients having diabetes mellitus or a weakened immune system (e.g., HIV positive, cancer chemotherapy, chronic steroid treatment, splenectomy)
- Patients who are bedridden (e.g., nursing home patient, stroke, chronic illness, recovering from surgery)
- A transplant patient (e.g., liver, heart, lung, kidney).

This thermometer is not intended for pre-term babies or small-for-gestational age babies.

This thermometer is not intended to interpret hypothermic temperatures.

Do not allow children to take their temperatures unattended.

Please consult your doctor if you see symptoms such as unexplained irritability, vomiting, diarrhea, dehydration, changes in appetite or activity, seizure, muscle pain, shivering, stiff neck, pain when urinating, etc., even in the absence of fever.

Even in the absence of fever, those who exhibit a normal temperature may still need to receive medical attention. People who are on antibiotics, analgesics, or antipyretics should not be assessed solely on temperature readings to determine the severity of their illness.

Do not modify this equipment without the authorization of the manufacturer.

Product description

- 1 Power button
- 2 Scanner
- 3 LCD display
- 4 Temperature button
- 5 Battery cover (2 X 1.5V type AA (LR6) batteries)
- 6 Protective cap
- 7 Silent mode switch
- 8 Guidance light

Why Braun No touch + forehead thermometer?

Measurement in under 2 seconds

The innovative infrared technology allows measurement of forehead temperature in under 2 seconds.

Accurate and reliable

Get the same professional accuracy whether touching your child's forehead or holding up to 5 cm away.

Easy to use

A measurement can be taken even while a child is sleeping.

Safe and hygienic

No-touch option helps minimize spreading of germs.

How does Braun No touch + forehead work?

The No touch + forehead thermometer measures infrared energy radiated from the skin at the center of the forehead area. This captured energy (which is twice as much thermal energy compared to a traditional forehead thermometer¹) is collected through the lens and converted to a body temperature value.

¹ A traditional forehead product without any optical system to capture radiated heat.

Obsah je uzamčen

Dokončete, prosím, proces objednávky.

Následně budete mít přístup k celému dokumentu.



Proč je dokument uzamčen? Nahněvat Vás rozhodně nechceme. Jsou k tomu dva hlavní důvody:

- 1) Vytvořit a udržovat obsáhlou databázi návodů stojí nejen spoustu úsilí a času, ale i finanční prostředky. Dělali byste to Vy zadarmo? Ne*. Zakoupením této služby obdržíte úplný návod a podpoříte provoz a rozvoj našich stránek. Třeba se Vám to bude ještě někdy hodit.

**) Možná zpočátku ano. Ale vězte, že dotovat to dlouhodobě nelze. A rozhodně na tom nezbohatneme.*

- 2) Pak jsou tady „roboti“, kteří se přiživují na naší práci a „vysávají“ výsledky našeho úsilí pro svůj prospěch. Tímto krokem se jim to snažíme překazit.

A pokud nemáte zájem, respektujeme to. Urgujte svého prodejce. A když neuspějete, rádi Vás uvidíme!