BRAUN

ThermoScan®

IRT 3030

GB	English	3
CZ	Česky	12
HR	Hrvatski	21
HU	Magyar	30
PL	Polski	39
RO	Română	49
RU	Русский	59
SI	Slovenščina	70
SK	Slovenčina	79
TR	Türkçe	88
SA	اللغة العربية	105

This product is manufactured by
Kaz Europe Sàrl under a license
to the 'Braun' trademark. 'Braun' is a
registered trademark of Braun GmbH,
Kronberg, Germany.

ThermoScan is a registered trademark
of Kaz Europe Sàrl

Kaz Europe Sàrl
Place Chauderon 18
CH-1003 Lausanne, Switzerland

© 2016, Kaz Europe Sàrl

www.hot-europe.com

Designed in Germany Made and printed in China


Patents: www.kaz.com/patents/braun

IRT3030EE

PN: 31IM303E190R1
25APR16

Printed in China

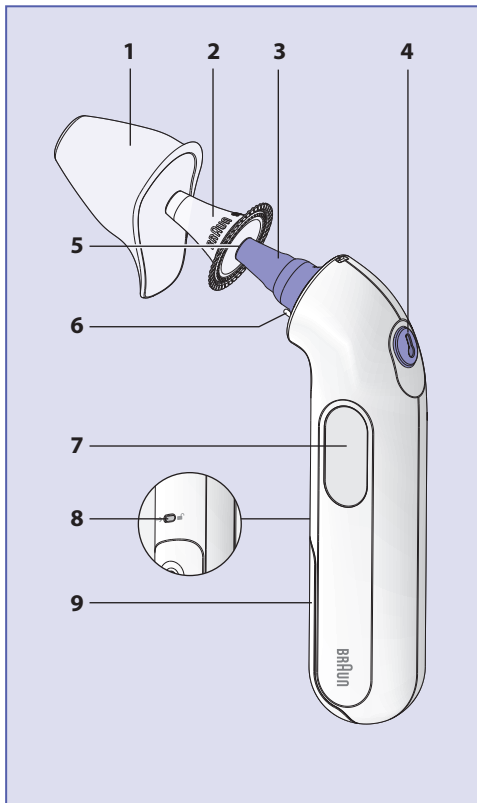Product Description

- ① Protective cap
- ② Disposable lens filter
- ③ Probe
- ④ «start» button
- ⑤ Probe tip (lens)
- ⑥ Lens filter detector
- ⑦ Display
- ⑧ Battery door safety lock
- ⑨ Battery door (1 X 3V type 2032 battery)

The Braun ThermoScan® thermometer is indicated for intermittent measurement and monitoring of human body temperature for people of all ages.

Important safety information

Please consult your doctor if you see symptoms such as unexplained irritability, vomiting, diarrhea, dehydration, changes in appetite or activity, seizure, muscle pain, shivering, stiff neck, pain when urinating, etc. in spite of absence of fever.

The operating ambient temperature range for this thermometer is 10 – 40 °C (50 – 104 °F).

Do not expose the thermometer to temperature extremes (below –25 °C / –13 °F or over 55 °C / 131 °F) or excessive humidity (> 95 % RH).

This thermometer must only be used with genuine Braun ThermoScan® disposable lens filters. Other lens filters can lead to inaccuracy.

To avoid inaccurate measurements always use this thermometer with a new, clean, disposable lens filter.

If the thermometer is accidentally used without a disposable lens filter, clean the lens after use (see «Care and cleaning»).

Keep disposable lens filters out of reach of children.

This thermometer is intended for household use only.

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