**TESLA SMART
COMPOSITION SCALE
WI-FI STYLE**
USER MANUAL

EN

DE

CZ

SK

PL

HU

SL

HR

UA

RU

Thank you for choosing our WI-FI scale. Our App allows users to easily connect their device to make health and weight control easy. Measurements are sent wirelessly to your mobile device and can be viewed using the **FREE app: Tesla Smart**.

Product Information:

- Max Capacity 180kg Division 100g
- Works with the “Tesla Smart” App and analyzes over 13 body health values.
Traditional values: Weight / BMI / Body fat (percentage) / water content (Body water percentage) / Muscle/Calorie (KCAL) / Bone Mass
New values: Body Score /Fat-free Body Weight / Visceral fat(grade) / Protein / Body Type / Metabolic Age
Curve trends lines for Weight / body fat / BMI
- App supports share to family multiple users
- Auto on & Off
- Low battery and overload indicator

Note: The scale will turn on when a weight of 5kg or more is placed on it. It is recommended to keep the scale in a working environment of: 0°C-40°C/32°F-104°F; Relative humidity<85% so that the health values will continuously be updated.

App “Tesla Smart” Information:

Search and download the free app “Tesla Smart” on Android or iOS.

System requirement for use of the App:

- iOS 8.0 or higher
- Android 4.4 or higher
- Due to the many different device manufacturers, models and software versions, full functionality cannot be guaranteed.

Please read the following instructions carefully before use.

The scales are designed for household use only, not to be used for commercial purposes.

Operation:

1. Open the battery compartment on the underside.
2. Insert batteries, please ensure that the polarity is correct.
3. Close the battery compartment again.

Note: If the message “LO” appears on the display, replace the batteries.

How to measure Body Weight

If you prefer to not connect to WI-FI and would just like to use the scale for simple weight measurement, then simply follow the steps below:


1. Press the platform center and remove your foot.
2. 0.0 will be displayed on the screen.
3. The scale will switch off and is now ready for use.

4. Step onto weight platform until your weight numbers appear on the screen.

Note: This initialization process must be repeated if the scale is moved. Always place the scale on a hard and flat surface to have accurate measurement. If you find a 3kg discrepancy 2 times consecutively, it is possible that the initialization process was not successful. Repeat steps 1-3.

Use with the App

To download the app, search for "Tesla Smart" in any app store. OR **scan below QR code**

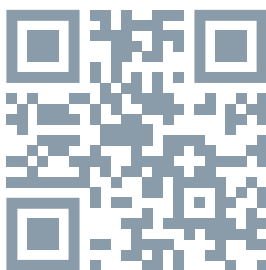Register your account:

follow the register page accordingly fill in the necessary information required.

Household / Member Management:

1, When logged in, Click right bottom  symbol, enter into "Home Management" to start "Add home"

2, Successful setting the home information, Click Top Left  symbol, to choose the one you set. Top Left corner will show the name.

Add Device:

1, In Welcome home page, click "Add Device" or click "+" symbol on the right top corner to enter into device pair page.

2, Choose "Exercise&Health" to find the device category "Body Fat Scale"

3, Long press the scale base housing Pair button till scale front display showed connection detecting square symbol and Wi-Fi signal symbol flashed fastly. Then click app "Confirm indicator rapidly blink" button.

Remark: if the WI-FI signal symbol flashed slowly, continue press the pair button until this symbol flashed fastly. Slowly flashed WI-FI symbol may failed WI-FI connection.

Pairing button

detecting square symbol and Wi-Fi symbol

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